



**VINNO**



**ZITRON MEDICOM SYSTEM (001321633-M)**

MEDICAL EQUIPMENTS SUPPLIES & SERVICES


Suite 3-6, Lower Level 3, The Horizon Annexe, Avenue 7, Bangsar South City,  
No. 8, Jalan Kerinchi, 59200 Kuala Lumpur.

TEL: +603 7499 3971 H/P: +6012 653 7892

Email: rameshitron@yahoo.com

**VINNO**

**VINNO A6 is a brand-new Color Doppler Ultra-sound System. Adhering to VINNO's simplified design style and extremely smooth work flow, it will bring you a distinct diagnostic experience.**

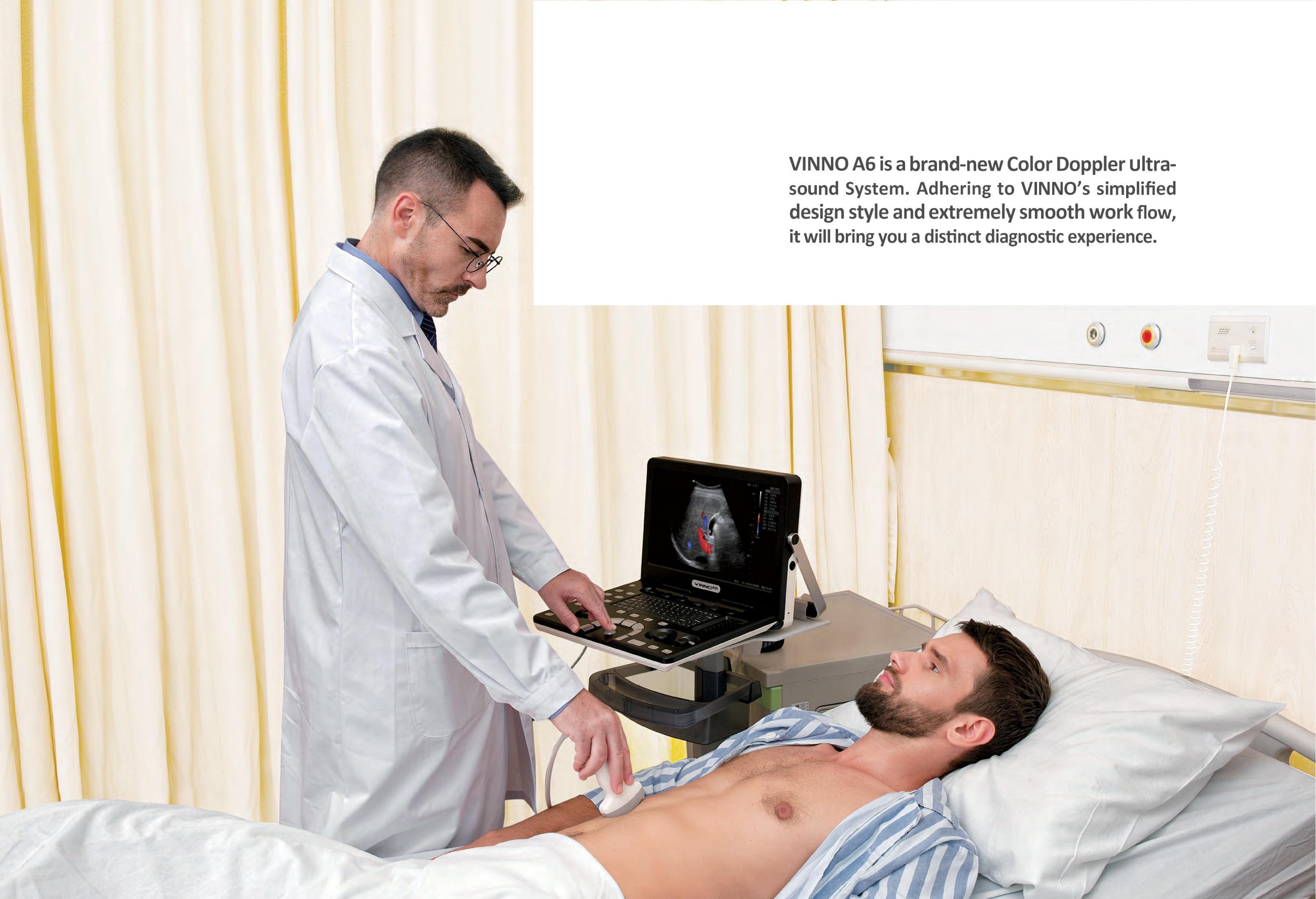


# Simplified Innovation

## 1. Simplified Platform

- » Minimize a powerful platform into a compact size to meet your daily work.
- » Innovative cooling solution delivers comfortable and silent working environments.



Intergrate

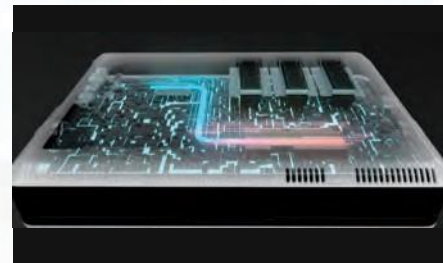

Compact

## 2. Simplified UI



**Fast Setting Area**-Complete common setting under scanning status.

**Outdoor and Indoor Solution**-Adjust monitor brightness according to the working environment.



**Quick Navigation Bar**-Browse long page fast by moving the trackball after click the interested area.

**Back ICON**-Return to the top of the long page with only one click.




**Simple UI**-Select preset directly without reiterative page turning.


**Magic Drag and Drop**-Adjust the sequence of measurement item easily by Magic Drag and Drop operation.

# Intelligent Workflow

## 1. Auto Image Optimization

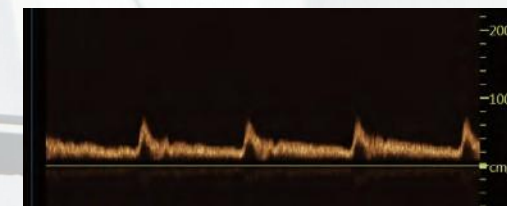

off

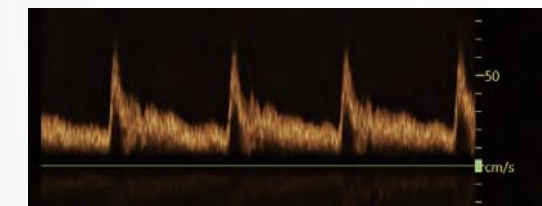

on

## 2. Auto Calculation


### Auto OB

Just one click, the system will recognize the target automatically and the gestational age will be displayed there.



### Auto IMT

Just one click, the system will display max, average values of intima thickness and valid points at the same time.

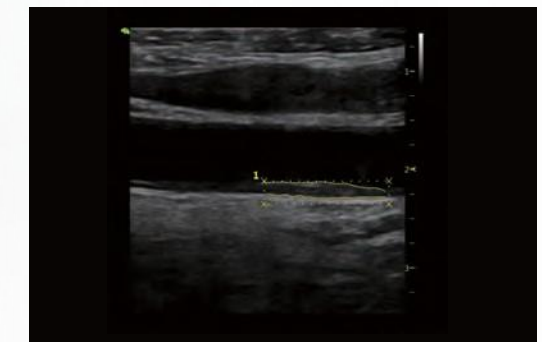




## Excellent Image Quality

**VFusion:** Compounding technology together with RF-platform, VFusion helps to improve the image information, enhances contrast resolution and reduces the noise by scanning from different angles and combining those images.

**VSpeckle:** Reduces speckles in homogeneous tissue and enhance edge sharper than conventional image processing.

**VTissue:** Matches different acoustic velocity with different tissue automatically. It ensures the acquired image with more uniformity and accuracy.

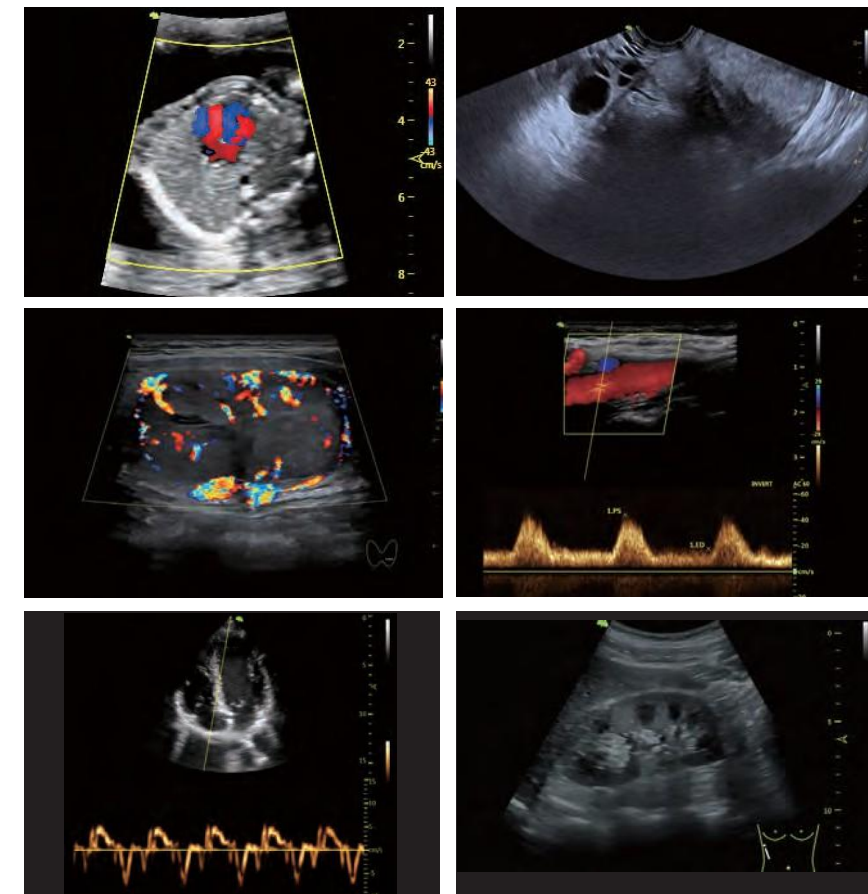

# Beyond Your Expectation

- » 15.6 inch LED display
- » Monitor can be tilted up to 40 degrees
- » Detachable probe holder
- » Support up to 3 probe ports
- » battery lasting for about 1.5 hours
- » Support CWD, ECG, 4D function
- » Support WiFi/Bluetooth connectivity
- » Support transducer: A2-5C, A4-12L, A1-4P, A4-9E, A3-6D





## Enhance Your Skill

"Education is the passport to the future". Log on VINNO's Flyinsono APP and VCloud platform to study continuously and share your experiences.

